



EquiRay

X RAY SYSTEM

Complete X-ray
system
for hippiatry
and large animals

Radiology on large animals or for equine care allows the veterinarian to provide a complete service to his clients.

The software contains a section dedicated to pre-purchase exams for horses.



- ① **100 kV - 30 mAs - 2kW**
Usable with standard socket
- ② **8,5 kg**
Lightweight for easy transport and use
- ③ **X-Ray Switch**
With 4 meters lenght
- ④ **Laser localizer**
Laser for correct centring
- ⑤ **Collimator**
43x43cm high brightness

The Porta 100 HF generator provides high reliability and with its power it is able to emit radiation doses sufficient to perform any type of examination.

Comes with an aluminum case for safe transportation and storage.
Extremely handy for manual use with the possibility of being mounted on a tripod.



- ① **Dimension 25x30 cm**
Standard dimension for hippiatry
- ② **Wi-Fi interface**
Fast connection to the computer
- ③ **Thin edges**
Useful for special exams on hoof
- ④ **OLED display**
Indicates the status of detector

The Vivix-V 2530VAW panel is a direct digital detector with wireless connection and cesium scintillator with 124 micron pixel pitch.


Long-lasting battery for up to 7.5 hours of uninterrupted use (8 hours in stand-by)
Fast connection to the software which guarantees a diagnostic image in about 2 seconds.



- ① **Laptop**
Lightweight powerful laptop
- ② **High brightness**
Clear vision in any situation
- ③ **VxVue software**
Full dicom for acquisition and post processing
- ④ **Wireless connection**
Fast connection to the detector
- ⑤ **Storage**
High capacity for studies storing

Software with the best post processing algorithms currently available on the market. The images returned by the digital system are processed and optimized in detail and in the balance of the black-white color curves.

The supplied hardware guarantees the best post-processing speed as well as a high storage capacity.

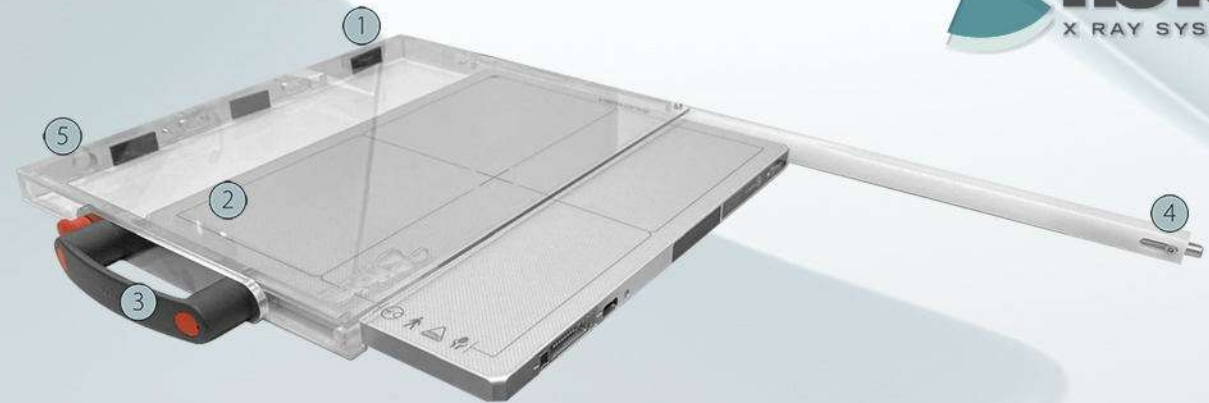

- ① **Detector protection**
Made in technopolymer for safety
- ② **Robust material**
For digital detector protection
- ③ **Repositionable handle**
Triple-sided for easy use
- ④ **Opening lock**
For more security
- ⑤ **Easy Access**
Easy access for charger and switch button

Technopolymer protection that guarantees both practicality and safety for the direct digital panel.

The repositionable handle allows you to rotate the panel ensuring a comfortable grip.

Options for skylines are also available for specific exams in horse medicine.

The digital system we offer includes all components needed for the execution of radiological exams.

EquiRay

X RAY SYSTEM



IBIS s.r.l. Via Cascina Bruciata, 3 Seriate, Bergamo
www.ibisray.it - info@ibisray.it
Tel.: 035 - 42 36 343