

Simply HP

X RAY SYSTEM



Simply is a lightweight and easy to handle mobile unit for maximum freedom of movement. Ideal for diagnostic activities in the ward and emergency room.

The low weight and simplicity of handling make it the most valid tool for hospital wards, emergency rooms and plaster rooms.

Radiological techniques

Two-point radiological technique with choice of kV and mAs or three-point radiological technique with choice of kV, mA and ms. All selected values are displayed on the touch screen display.

Anatomical Techniques


It is possible to program and recall up to 240 tests via the virtual keyboard on the touch screen.




Dosimetric system

Simply can be equipped with DAP. Patient data can be entered via the virtual keyboard; dose data can be printed for convenient archiving.

Automatic shutdown

Simply is equipped with an automatic shutdown system that turns off the equipment after 30 minutes of inactivity to prevent damage to the insert.

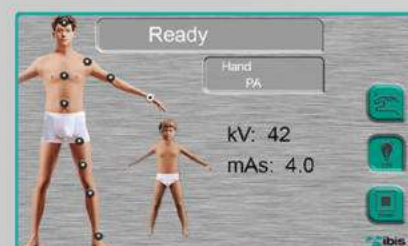

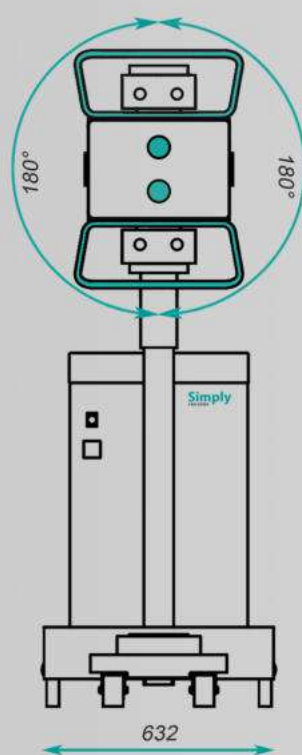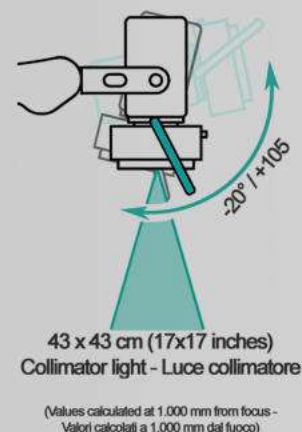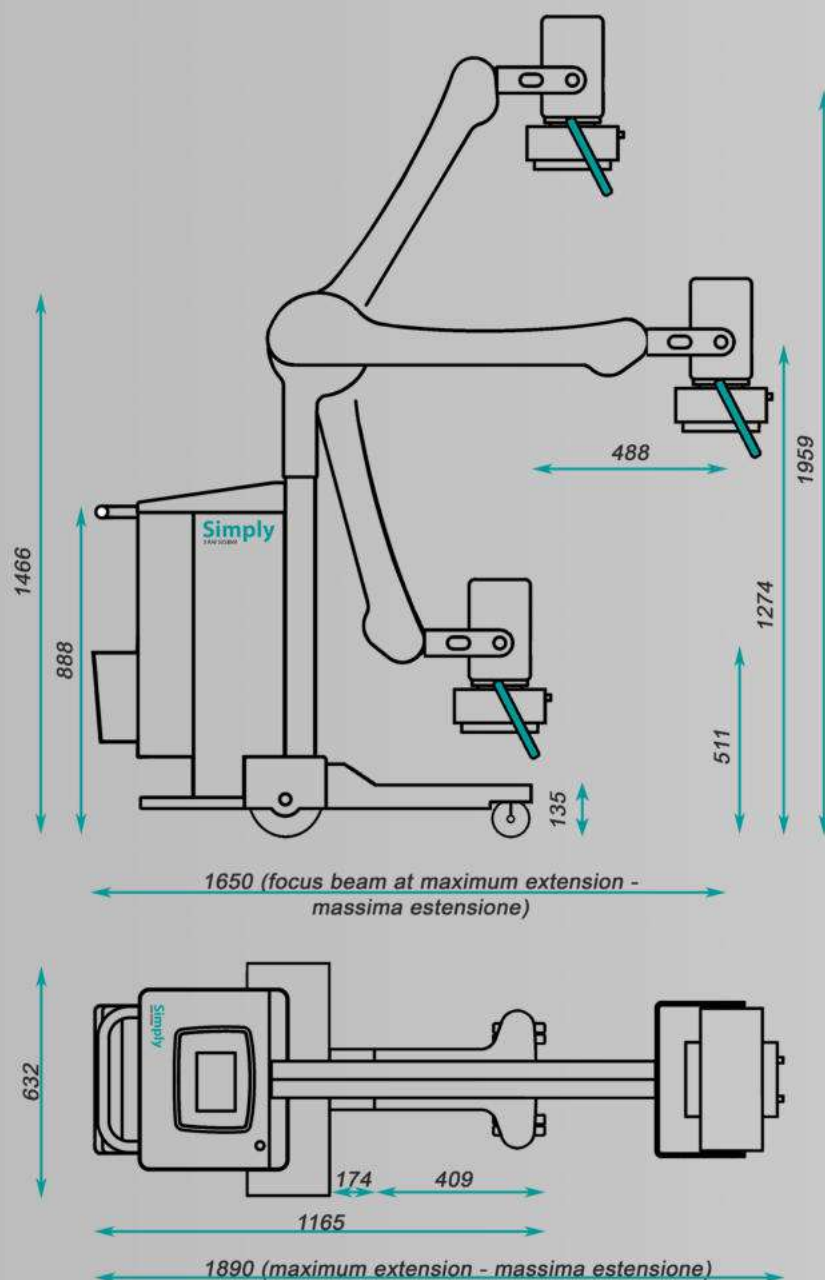

Technical features

Power:	15 kW
kV range:	40-125 kV
mAs range:	0,2-220 mAs (28 steps)
Focal Spot:	0,6-1,3 mm
Thermal capacity:	900 kJ - 1024 KHU
Max mA:	320 mA
Frequency:	100 kHz
Power supply:	230Vac $\pm 10\%$ 50/60 Hz

Simply HP15

Simply HP32

Power:	32 kW
kV range:	40-125 kV
mAs range:	0,2-220 mAs (28 steps)
Focal Spot:	0,6-1,3 mm
Thermal capacity:	900 kJ - 1024 KHU
Max mA:	400 mA
Frequency:	100 kHz
Power supply:	230Vac $\pm 10\%$ 50/60 Hz

IBIS s.r.l. has been designing and producing medical equipment for radiology since 2001, maintaining high quality standards throughout the entire production process.

The company specializes in the production of mobile ward units, "C"-arms and veterinary radiological equipment.

The product range is completed with the distribution of portable x-ray systems for home use and with the range of digital detectors.



CE 0051



IBIS s.r.l.

Headquarter: Via cascina bruciata, 3 24068 Seriate, Bergamo - Italy
www.ibisray.it - info@ibisray.it - Tel.: +39 (0)35 4236343

www.ibisray.it



Cod. 3582